


1.5 mL microfuge Tube/ 1 mL Tube/ PCR Strip Magnetic Separation Rack


**All Permagen Labware
Magnetic Separation Racks are
designed for both performance and
durability.**

- Separate magnetic beads from up to (24) 1.0 mL screw cap tubes which are spaced apart to match any 96- well microplate or storage box allowing simple, direct transfer
- Separate in one or two PCR strips
- On reverse side, separate in up to (12) 1.5 mL microfuge tubes. Holes are angled directing tubes to come in close proximity with the magnets
- Strong Neodymium magnets for fast separations of many magnetic bead chemistries
- Solid, one piece aluminum construction eliminating the concern for material inconsistencies or breakage associated with some plastic products on the market for trouble free use
- Hard-coat finished to maximize durability
- Includes two clips which will aid in pushing PCR strips close to magnets (see last image)
- Volume discounts available
- Most durable separation racks on the market at competitive prices
- *Made in USA*

For more information please visit our website or email us at info@permagenlabware.com

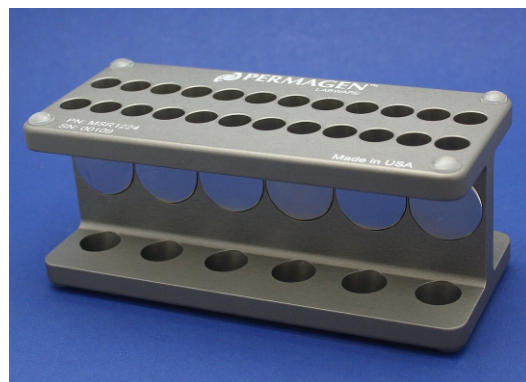

MSR1224 Magnetic Separation Rack

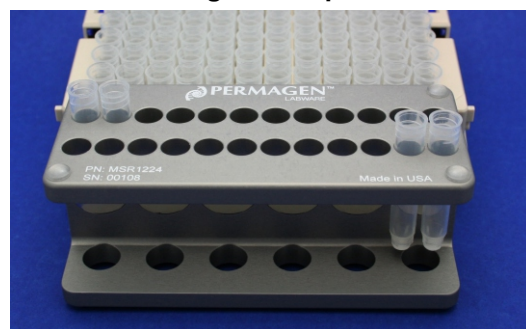

Simple transfers between 96- Well



One rack, several options

