


## U Series 96- Well Magnet Plate


All Permagen Labware Magnetic Separation Plates are designed for both performance and durability.

- Solid aluminum construction for precise, trouble free use eliminating the concern for material inconsistencies or breakage associated with some plastic products on the market
- Precision engineered, giving you repeatable, reliable results
- Unique, fastener free construction, equals no screws or pins to come loose
- Side clearance for liquid handler grippers
- Hard-coat finished to maximize durability
- 96 Strong Neodymium ring magnets for fast separations of many magnetic bead chemistries  
\*beads are pulled into a small ring (approx. 4 mm) near well bottoms
- Integrated raised corners, making it compatible with almost any Microplates and PCR plates.
- SBS compliant footprint top and bottom (raised bottom to accommodate platforms with deeper nests)
- Manual or robotic use
- Available in three magnet strengths to suit your needs  
U500 – Up to 2 mL max or smaller volumes faster  
U480 – Up to 1 mL max or smaller volumes faster  
U380 – Up to 350 µl
- Volume discounts available
- *Made in USA*

For more information please visit our website or email us at [info@permagenlabware.com](mailto:info@permagenlabware.com)

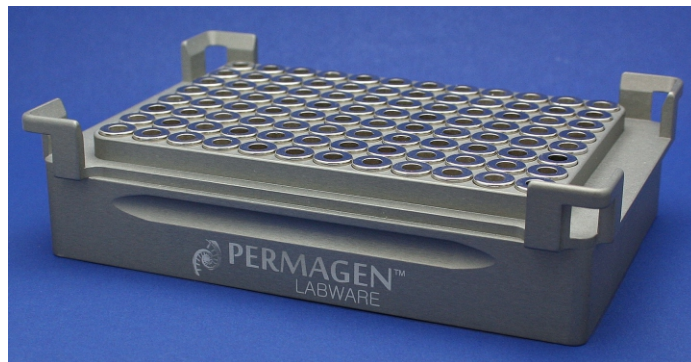

**U series 96- Well Magnetic Separation Plate**

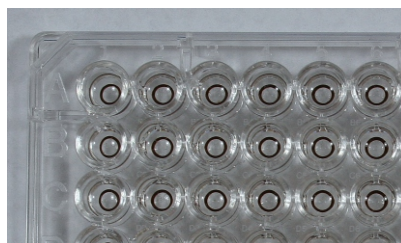

***Magnetic Beads are pulled into rings near the bottom of each well***



**Compatible with almost any 96- Well Microplates and many liquid handling platforms**